


# Rechargeable Sonic Toothbrush Instruction

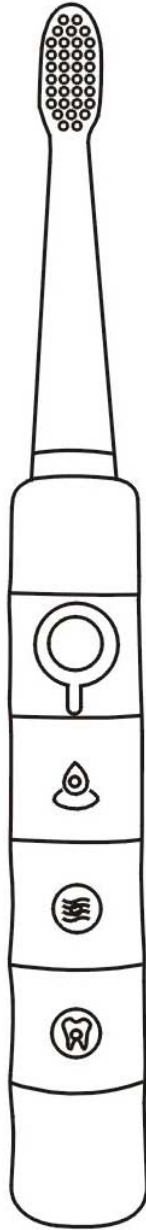

Before operating this unit, please read this

Thanks for purchasing our rechargeable sonic toothbrush very much

This operating instructions contain the safe operation and instructions for proper use of this product. Please read this manual carefully before use and be careful to keep it from damage and danger

### Safe operating instructions

1. Please do not use it outside the specified use
2. When connected with power, please check whether the power voltage is suit for your local voltage or not,
3. Please use original charger to avoid dangerous
4. Please avoid abrasives with whitening and removing tooth stains to use together
5. Please do not high strength use for teeth and gums
6. Children and can not operate their own or the mouth feel weaker, please do not use
7. If the teeth are being treated, or if the implanted tooth is not strong, do not use it
8. Do not drop it, do not hit it violently
9. Do not soak in water for long periods
10. This product is only for personal use and is not available to dental clinics or organizations for multiple patients


### warn

Children and persons with physical and mental disabilities are required to use under the supervision of the guardian

In the case of periodontal disease and oral symptoms, please use under the guidance of a doctor

If you feel uncomfortable or hurt in your mouth, please stop using and consult dentist.

When the electric toothbrush fails, please consult your local government how to dispose of this tool to maximize the protection of the environment

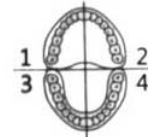

Four oral areas


### How to use

1. Wet the toothbrush head, and squeeze appropriate toothpaste on it.
2. Place the bristles on the teeth of the 45 degree angle aligned with the gum line.
3. Press ON/OFF button, brush start working or stop working.
4. You just gently move the bristles to place where need to clean, The pressure should be within 150g (Sonic toothbrush is to rely on high-frequency vibration of the bristles to clean the mouth, so when brushing teeth, the using force teeth is much smaller than ordinary, about 1/8 or so can)
5. Dentist advice that the oral cavity is divided into 4 areas(as shown above), slowly brushing the toothbrush during the brushing process to ensure that each area has about 30 seconds of cleaning time, the machine has a 30-second area reminder function, every 30 Seconds will be suspended for 0.5 seconds to remind you to replace the area to brush your teeth.,
6. Press on/off button, brush start working, then press button keeping three seconds to stop working. This toothbrush has six brushing modes: clean (Strong/weak), soft(Strong/weak), whiten(Strong/weak). You can press the one to work mode selection to own needs. This toothbrush also has a memory function working mode, next time you use the brush, directly into the mode you has used, no need to adjust.



## Product Features

1. High-frequency sonic vibration: 31,000 times per minute high-frequency acoustic vibration, just 15 days, you can remove the stains (including tea stains, smoke stains), to achieve whitening effect
2. Ergonomic design: waterproof non-slip, feel comfortable, safe and reliable
3. USB fast charge design ,Unplug the waterproof plug at the bottom of the toothbrush. Plug the DC plug of the charging cable into the charging port at the bottom of the toothbrush before inserting it, and then insert the USB plug into the computer or the USB plug on the 5V. The maximum charging time is 6 hours. The green light flashes lowly when charging, and the green light is on when charging is full.
4. Three brush head: standard with three soft brush heads, can be replaced anytime, anywhere, long-term care of your oral health
5. Safety brush cap more healthy: configure the security brush cap, not only to prevent bacterial invasion, more conducive to carry convenient and healthy

## Name of each part

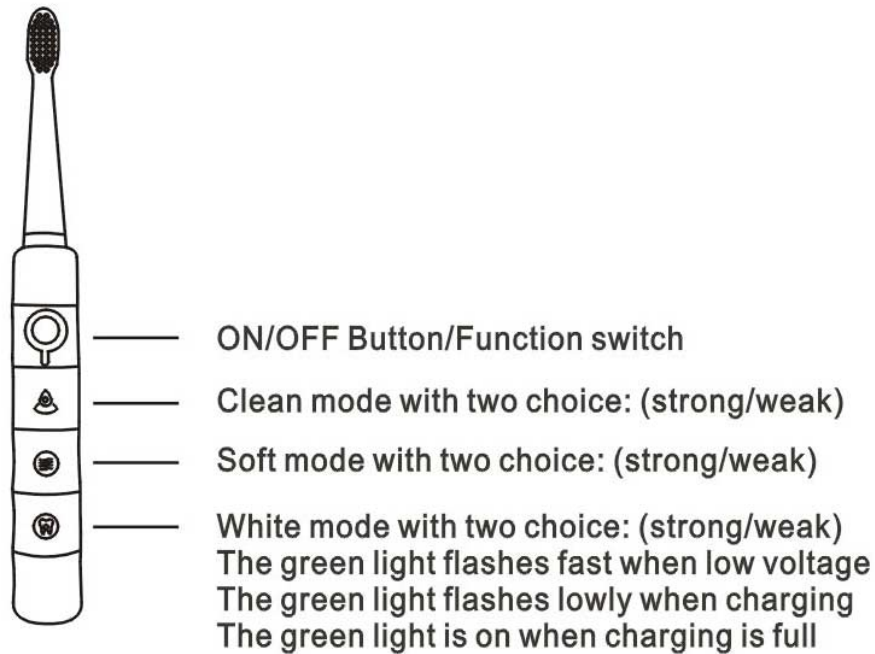

## Brushhead installation

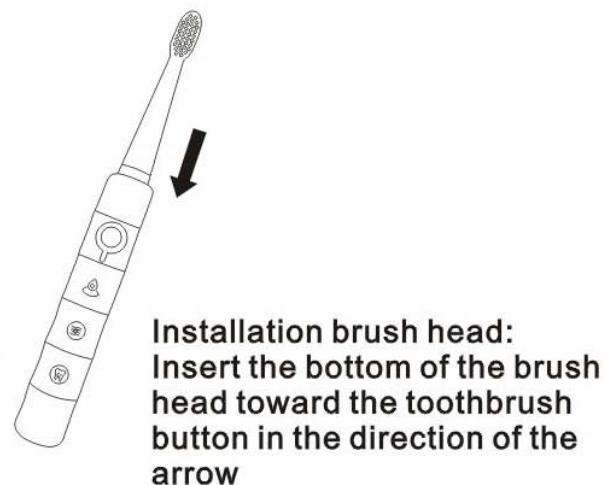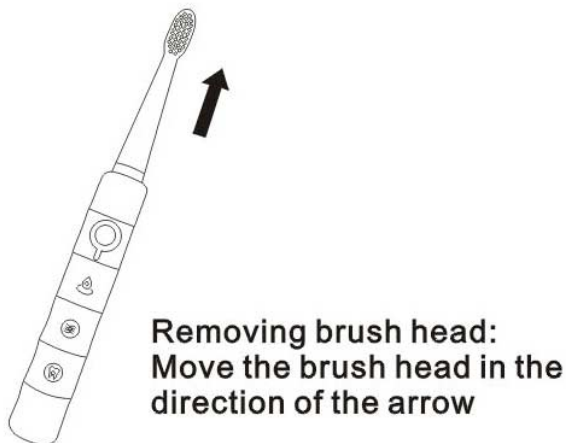

## Specifications and implementation standards

Parameter	specification
Battery is full (The indicator light is always on)	Weight: Toothbrush handle about 56G
Rechargeable battery: lithium-ion 500Ma	Size of toothbrush: 208*24*24MM
Rated voltage: DC5V 50/60Hz	Six modes: Clean(strong/weak), Soft(strong/weak), White(strong/weak)
Rate current: 3.7V 300Ma	Bristles: Wavy DuPont soft Nylon Bristle
Motor speed: 16500 rpm / min	High-frequency sonic vibration: 31000 times/min
Rated power: 1W when charging	

## Accessory



Soft toothbrush heads 3



Dustproof cover 1



USB charging cable 1

## Environmental protection

The outer packaging consists of environmentally friendly raw materials, so it can be discarded in the general trash can.

When the electric toothbrush is aged, consult the local government how to handle the tool to maximize the protection of the environment